


AUTO CHART PROJECTOR

ACP-8_RNEW

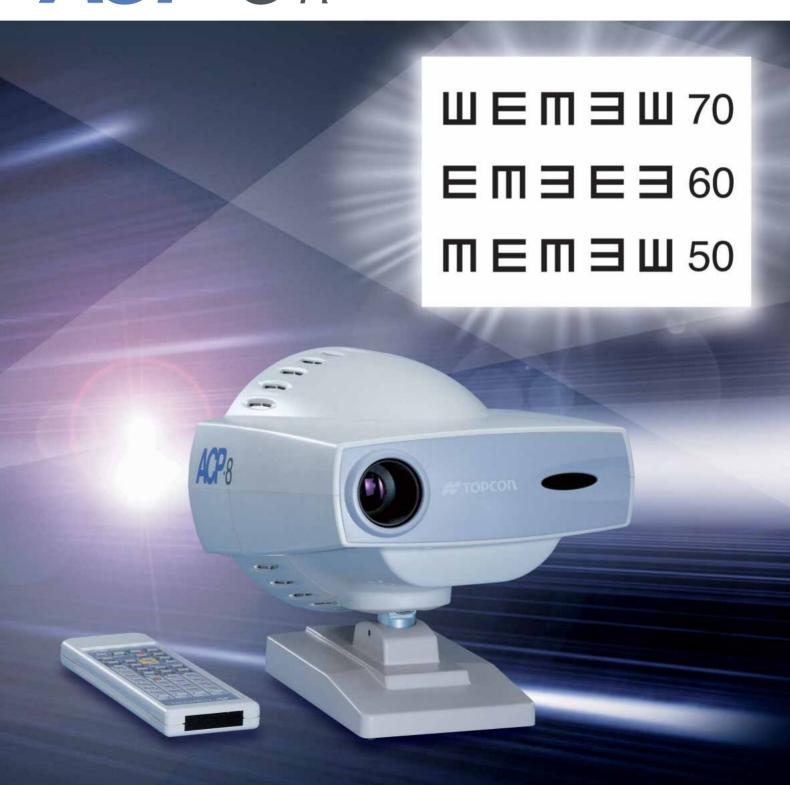

VERSATILE...

QUICK... CONVENIENT...
FUNCTIONAL

The TOPCON ACP-8 R Auto Chart Projector provides improved sharp and clear image projection and fast chart rotation.

A cordless remote controller allows

A cordless remote controller allows instant projection of the 30 AO compatible test charts at a rate of 0.03 seconds per frame. In addition, the ACP-8 $_{\rm R}$ can offer 3 $_{\rm X}$ 5 character charts with a wide projection size (330 $_{\rm X}$ 270) for fast and efficient testing.

Features

Two Program Settings

The ACP-8 R can be programmed to project a specific sequence of charts for maximum flexibility and refraction ease.

Two programs can be memorized.

Single Character Masking

The versatile masking capabilities of the ACP-8 $\,{}_{\rm R}$ allow projection of a single character at a time.



Easy Setup and Maintenance

Focusing setup is easily adjusted without opening the cover. It is also very easy to replace the bulb. You can replace the bulb without any adjustment.

Coupled with CV-3000

The ACP-8 R can be remote controlled directly with the TOPCON Compu Vision CV-3000 remote controller. This combination allows the operator to quickly and easily access all functions for general visual examination.



29310

856

Tilt Function

The horizontal axis of the projector can be adjusted to guarantee a perfectly positioned projection on the screen.

Variable Focus Lens

The variable focus lens of the ACP- $8_{\,\rm R}$ allows the instrument to be positioned wherever it is most convenient.



H Type chart



ENH 200 **SLC** 150

H B V 100 P H T 80 V L N E A 70 D A O F C 60 E G N D H 50 F Z B D E 40 O F L C T 30 A P E O F 25 TZVEC 20 OHPNT 15 EVOTS 10

EVOTL 20 TBGAB 20 HNFZC 20 APLBK 25 6 3 8 2 5 20 ZRHCT 15

ПЭШ 100 ЭШЕ 80 ШЕМ Э Ш 70 ЕМ Э Е Э 60 М ЕМ Э Ш 50 Ш Э М Е Ш 40 Э М Е Ш Э 30 М Е Ш Э М 25 E Ш 3 M 3 20 M Ш E 3 M 20 Ш 3 M E Ш 15

100

 발생하를 **~** 40 라 를 발하성 30 참**~** 선플라 20 9 3 3

9375840 3829430 5394625



PT TP ZN NZ DA AD TH HT

587₂₀₀
426₁₅₀

293₁₀₀ 856₈₀ 3845270 9683760 4257650 8257320 4783615 5864910













Technical Specifications

Refracting distance: 2.9 to 6.1 m **Variable focus projection distance:** 2.9 to 7.1 m

Projection magnification:30x (in 5m refraction)Projection size: $330 \times 270 \text{mm}$, $\phi 300 \text{mm}$

Number of charts: 30

Chart change-over: 1 frame/0.03 sec.

Number of masks: Open 1, Horizontal line 5, Vertical line 8, Single isolation 21, R & G 1

Mask change-over: 1 frame/0.02 sec.

Program step: 2 type, Max. 30 steps available
Projection lamp: 12v 50w (Halogen lamp)

Automatic shut-off: After 10 minutes

Electricity: AC 120, 220, 230 or 240V, 50/60Hz

Power consumption: 80VA

Dimensions: 226 (W) x 300 (D) x 245 (H)mm *including Base height

Weight: 6.0kgs

Standard Accessories

Polarized Metal Screen

(21.7" x 17.7" with non-glare finish)



Optional Accessories







*Subject to change in design and/or specifications without advanced notice.

IMPORTANT

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

-